VIR-ELISA ANTI - RUBELLA- IgM (Immunocapture assay)

1. INTENDED USE

Immunocapture assay for qualitative determination of IgM specific antibodies to Rubella virus in human serum and plasma.

IgM REF EM 126 ▽ 96 IVD

2. PRINCIPLE OF THE ASSAY

The ELISA test ANTI - RUBELLA - IgM is based on the anti- μ -capture principle. Anti-human IgM antibodies are fixed to each well of the microtiterstrips. All the IgM class antibodies present in the patient's sample are bound during the first incubation to the solid-phase antibody. After removing unbound material by washing rubella antigen/POD conjugate is added. During the second incubation the conjugate will bind to the Rubella specific IgM antibodies, which have been captured by the anti-human IgM antibodies during the first incubation. Excess conjugate is then removed by washing and TMB substrate is added, resulting in the development of a blue colour. The enzyme reaction is terminated by the addition of a stop solution. The intensity of the yellow colour thus developed is proportional to the concentration of antibodies in the sample.

3. DIAGNOSTIC RELEVANCE AND INTERPRETATION OF RESULTS

Shortly after the viraemic phase (which lasts approx. 5-6 days) antibodies are generated. Acute infection results in high IgG as well as IgM antibodies. The detection of both antibodies in high concentration is possible for up to 8 weeks, then the level of antibodies falls. While IgM antibodies disappear generally, IgG antibodies persist life-long.

Following a Rubella immunisation the increase in antibody is significant, although the concentration of antibodies is generally lower than for an infection. Specific IgG antibodies can be detected for up to 15 years, IgM on average 8-12 weeks. In special cases, persisting IgM antibodies were found over several months.

Serological findings in a Rubella virus infection are strongly dependent on the stage and duration of the clinical symptoms. A single antibody activity, even a high level activity, does not provide proof of a recent infection and cannot necessarily be equated with adequate immune protection. To obtain a final diagnosis the patient history and clinical symptoms should be taken into consideration.

Note!

In case of a positive IgM antibody detection it is essential to confirm this result with at least two other different tests !

IgG	IgM	Interpretation	Recommendation
-	-	no specific antibodies (Ab) detectable, however infection possible	By suspect of acute infection – control tests are recommendable (after 7 days at the latest)
+	-	probable previous infection or vaccination, reinfection possible	Monitoring of IgG antibodies (Sample collection within 10-14 days), a significant increase of IgG Ab concentration in the absence of IgM Ab indicates a reinfection
-	+	primary infection probable	Monitoring of IgG and IgM Ab to determine seroconversion; supplementary tests e.g. IFA,CF,HAH and determination of the avidity of the IgG Ab.
+	+	acute infection probable, vaccination possible	Monitoring of IgG and IgM antibodies; supplementary tests e.g. IFA,CF,HAH and determination of the avidity of the IgG Ab

- negative, + positive

4. PERFORMANCE CHARACTERISTICS

Specificity / Sensitivity

292 samples were tested in parallel in VIR-ELISA ANTI-RUBELLA-IgM (Immunocapture assay) and comparison methods (ELISA). The specificity and sensitivity are based on the results found.

Specificity: IgM 99,2% **Sensitivity:** IgM 100%

For the calculation of the **specificity and sensitivity** of the ELISA test, equivocal results were defined as positive. The results refer to the groups of samples investigated.

Precision and reproducibility

Intra-assay reproducibility was determined by testing samples of different levels of antibody reactivity for at least 22 times in one test run. The coefficient of variation (CV) of the reactive IgM samples was < 10%.

Inter-assay reproducibility was determined by testing samples of different levels of antibody reactivity in 10 different test runs. The CV of the reactive IgM samples was < 10%.

Cross reactivity

A patient may have developed specific rubella IgM antibodies as a result of a polyclonal response to a primary infection such as with Epstein-Barr virus (EBV), Cytomegalovirus (CMV), Hepatitis A or even Gram-negative bacteria.

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6. KIT COMPONENTS

- MICROTITERSTRIPS** MTS
One microtiterplate is supplied, which contains 12 microtiterstrips of 8 breakapart wells. Each well is coated with mouse antibodies to human IgM.
- SAMPLE DILUENT** SPE|DIL|RUBELLA|IgM
One bottle containing 100 ml (or 2x50 ml) of ready-to-use sample diluent buffer Rubella IgM. The buffer includes 0.05% microbial preservative.
- PEROXIDASE CONJUGATE** CONJ|POD
One vial containing 12 ml of ready-to-use rubella antigen/POD conjugate with 0.049% Thimerosal as preservative.
- NEGATIVE CONTROL** CONTROL|-
One vial of 1.2 ml containing human serum with 0,095% sodium azide as preservative. Ready to use.
- CUT OFF CONTROL** CUTOFF
One vial of 1.2 ml containing human serum with 0,095% sodium azide as preservative. Ready to use.
- POSITIVE CONTROL** CONTROL|+
One vial of 1.2 ml containing human serum with 0,095% sodium azide as preservative. Ready to use.
- TMB SUBSTRATE** SUBS|TMB
One vial containing 13 ml of ready-to-use tetra-methylbenzidine (TMB) substrate.
- WASH SOLUTION 25X** WASH|BUF|25x
One bottle containing 80 ml (or 2x40 ml) of wash solution concentrate.
- STOP SOLUTION** SOLN|STOP
One bottle containing 15 ml of 0.95N H₂SO₄ stop solution. Ready to use.
- ADHESIVE SEALS**

The safety data sheet (MSDS) is available upon request.

7. STORAGE AND STABILITY

Store all reagents at 2-8°C. Protect them from intense light and do not freeze. The expiration date of each component is indicated on the respective vial label. Do not use reagents beyond the expiration date.

After opening, MTS must be stored at 2-8°C in the plastic bag with desiccant and are stable up to 4 weeks.

The diluted WASHBUF is stable up to 4 weeks if stored at 2-8°C.

Use only MTS with an intact vacuum packaging.

8.

MATERIALS AND EQUIPMENT REQUIRED BUT NOT PROVIDED

- Test tubes for sample dilution
- Timer
- Micropipettes, multipipettes 10-1000 μ l
- One-litre graduated cylinder, distilled water
- ELISA washer or multichannel pipette
- Spectrophotometer for micro-plates (450-nm/ reference wavelength 630/620 nm)
- Paper towels, pipette tips
- Incubator for \bar{x} 37°C with high relative humidity

9. WARNINGS OR PRECAUTIONS SAFETY PRECAUTIONS

The ELISA test is for **IVD** use only.

- Do not mix lot specific reagents, such as **MTS**, controls and **CONJ|POD** from different kit lots. **SUBS|TMB** and **SPE|DIL|RUBELLA|IgM** don't have to be from the original test kit. But the lot of these reagents should be the same as indicated on the kit label. **WASHBUF|25x** and **SOLN|STOP** can be used for all ELISA tests.
- Seal all bottles properly after use in order to avoid bacterial contamination. All samples and kit components should be considered potentially infectious. All controls have been tested for Hepatitis Bs antigen, anti-HIV I, anti-HIV II, anti-HCV (CE/FDA), and found to be negative.
- Do not pipette by mouth.
- Avoid contact with skin and mucous membranes when handling reagents, which contain preservatives (see kit contents). Wash thoroughly with water in case of contact and possibly look up a doctor.
- Controls containing sodium azide may react with lead and copper plumbing, building up explosive metal acids. Flush with sufficient water when disposing of reagents.
- The **SOLN|STOP** (0,95 N H₂SO₄) contains sodium hydroxide which may irritate skin and mucous membranes. Wash thoroughly with water in case of contact.
- For disposal the legal regulations have to be followed.

10. SPECIMEN COLLECTION AND STORAGE

- Only qualified and well-trained employees should carry out the assay procedure.
- The instruction for use describes the applicable test method. In case of modification or applications others than the intended use, or the use of automatic processors, the user has to validate the procedure and take the responsibility for it.
- Microbial contaminated specimen may cause interferences.
- Lipaemic, hemolytic or icteric samples should only be tested with reservations although in our testing no negative influence has been found.
- Suitable specimens are serum or plasma (heparinized citrated, EDTA) samples obtained by standard laboratory techniques.
- The samples should not be heat-inactivated since non-specific results may occur.
- Results on tests using CSF are not available.
- Patient samples should be stored at \pm 2-8°C. For long-term storage \pm -20°C or lower is recommended. Avoid repeated freeze-thaw cycles.
- Note:** Diluted patient samples must be used on the same day.

11. ASSAY PROCEDURE REAGENT PREPARATION

Bring all reagents to room temperature prior to use!

WASHBUF: Dilute the **WASHBUF|25x** 1:25 with distilled water e.g. add 40ml of **WASHBUF|25x** to 960 ml distilled water and mix well.

Dilution of samples: Dilute patient samples **1:101** with **SPE|DIL|RUBELLA|IgM** e.g. 10µl sample + 1ml **SPE|DIL|RUBELLA|IgM**, mix thoroughly.

Controls are ready to use !

Note: To take into consideration pipetting time, it is recommended to repeat the **CUTOFF** after every 6 **MTS** (resp. after a pipetting time of \geq 7 min.) to evaluate the following patient's tests with the new calculated cut-off value.

Take the required **MTS** out of the foil packets and place them in the holder. Possibly remaining wells of a **MTS** have to be stored at 2-8°C tightly sealed in the plastic bag provided, with the desiccant inside.

12. PIPETTING AND INCUBATION STEPS

- Pipette 100µl of each control and each diluted patient sample into the wells. **Leave a well (A1) empty for the blank (100µl SUBS|TMB).**
- Cover the **MTS** with an adhesive seal and incubate the wells for 60 minutes at 37°C (\pm 1°C) in an incubator with relative high humidity or in a humid chamber.
- Wash the wells four times as described in section k. **WASHING PROCEDURE**
- Add 100µl of ready-to-use peroxidase conjugate to each well **except of the blank (A1).**
- Cover the **MTS** with an adhesive seal and incubate the wells for 30 minutes at 37°C (\pm 1°C) in an incubator with relative high humidity or in a humid chamber.
- Repeat washing as in section C above.
- Add 100µl of ready-to-use TMB substrate to each well.
- Incubate the wells at room temperature (18...25°C), in the dark for 30 minutes
- Add 100µl of stop solution to each well. Tap gently to ensure homogenous colour distribution and read within 10 minutes.
- To read the plate, make sure that the bottom is free from moisture and that no air bubbles are in the wells. Read the absorbance of the well contents at 450nm on a suitable plate reader. On readers equipped with a dual wavelength facility set the reference filter to 620/630 nm.

Attention: The absorbance (OD) of the Blank must always subtracted from the OD values of the controls and samples.

PROCEDURAL NOTES

Do not allow the wells to dry out between incubations. Comply with the given incubation temperatures and times.

k. **WASHING PROCEDURE**

The washing procedure can be done manually with a multichannel pipette or on an automatic plate washer. Empty the wells, invert and tap dry on paper towel. Wash four times with a soaking time of approx. 30 seconds (300 µl).

13. SUGGESTIONS FOR TROUBLESHOOTING

In case that the ELISA instructions are followed strictly, the reagents are handled with care and the samples and reagents are pipetted carefully, the following kinds of errors can be avoided to a large extend.

ERROR	POSSIBLE CAUSES
No colorimetric reaction after addition of Peroxidase conjugate and TMB substrate	No Peroxidase conjugate dispensed, contamination of Peroxidase conjugate (possibly with control sera during pipetting) may cause an inactivation.
Generally too high reaction	Incorrect Peroxidase conjugate (i.e. not from original test kit), incubation time too long or incubation temperature too high, water quality for Washing Solution insufficient (low grade of deionisation).
Generally too weak reaction	Incorrect Peroxidase conjugate (i.e. not from original test kit), incubation time too short, incubation temperature too low.
Reagent blank too high	Incorrect pipetting of sample diluent, contaminated reagents, reagents expired, exceeding of incubation time and temperature, external contamination of the bottom of microtiterstrips, (clean carefully!).
False positive / negative samples	Incorrect dilution of samples, microbial contaminated specimen
Unexplainable outliers	Contamination of pipettes, tips or containers or with metals (iron, copper etc.), insufficient washing.
High variation (within a series)	Reagents (including microtiterstrips) not pre-warmed to room temperature prior to use. Washer is not washing correctly!
High variation (from series to series)	Incubation conditions not constant (time, temperature) high variation of incubation temperature, controls and samples are not carried out at same time (same intervals) check pipetting order, person related variation, strips dried out after washing (results are not reproducible).

14. VALIDITY OF THE ASSAY

All controls should be carried out with every test run.

The test must comply with the following validation criteria:

- OD-value of the Negative Control should be < 0.150,
- OD-value of the Cut-off Control should be >0.200,
- Ratio of the Positive Control/ Cut-off value should be ≥ 2.0 ,
- OD-value of the Blank should not be higher than 0.100.

If controls give invalid levels then results from the test samples are invalid too and retesting is required.

15. CALCULATION OF RESULTS

A QUALITATIVE CALCULATION

Calculation of "Cut-off Value"

CUT-OFF VALUE = OD VALUE OF THE **CUTOFF**

CUT-OFF RANGE = CUT-OFF VALUE \pm 10%

Interpretation of sample results:

RESULT	DEFINITION
negative -	OD value sample < Cut-off value -10%
equivocal	OD value sample \geq Cut-off value -10% OD value sample \leq Cut-off value +10%
positive +	OD value sample > Cut-off value +10%

Equivocal results should be retested. Following the confirmation of the equivocal result the monitoring of the patient's antibodies is recommended in order to exclude unspecific reactions, which may also cause equivocal results.

B CALCULATION OF RATIO (CUTOFF INDEX, COI):

Patient samples may also be quantified and interpreted by means of the calculation of the ratio (Cutoff Index, COI):

COI = OD value of sample/ Cut-off value,

whereby a ratio of 1.000 is equivalent to the Cut-off value.

Interpretation of sample results:

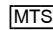
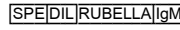
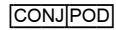
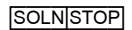
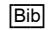
Ratio > 1,1 positive result
Ratio < 0,9 negative result
Ratio 0,9-1,1 equivocal result

Regarding diagnostic relevance and interpretation of results see **page 1**.

For further information please visit our website:

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