ELISA ANTI - VZV IgG/ IgM/ IgA

1. INTENDED USE

Enzyme-immunoassay for determination of antibodies against Varizella Zoster Virus in human serum and plasma.

IgG [REF] VEG 106 s 96 [IVD]
 IgM [REF] VEM 106 s 96 [IVD]
 IgA [REF] VEA 106 s 96 [IVD]

2. PRINCIPLE OF THE ASSAY

The detection of antibodies is based on the principle of an enzyme-linked immunosorbent assay (ELISA). The purified, homogeneous antigen is fixed to each well of the microtiterstrips. Any specific antibodies present in the patient's sample are bound during the first incubation. After removing unbound material by washing, the presence of specific antibodies is detected using Anti-human IgG/ IgM/ IgA conjugate during the second incubation. Excess peroxidase conjugate is then removed and TMB substrate is added, resulting in the development of a blue colour. The enzyme reaction is terminated by the addition of a stop solution. The intensity of the yellow colour thus developed is proportional to the concentration of antibodies in the sample.

3. DIAGNOSTIC RELEVANCE AND INTERPRETATION OF RESULTS

To obtain a final diagnosis the patient history and clinical symptoms should be included for the interpretation of the serological results and possible cross-reactivity should be taken into consideration.

| IgG | IgM | IgA | Interpretation | Recommendation |
|-----|-----|-----|--|---|
| - | - | - | No specific antibodies detectable. | By suspect of acute infection . control tests are recommended. |
| + | + | + | Acute or recent infection prob., reactivation or recent immunisation possible. | Monitoring of IgG, IgM and IgA antibodies (Sample collection within 10-14 days): confirmatory tests e.g. IFA, CF. |
| + | - | + | | |
| + | + | - | | |
| + | - | - | Previous infection prob., reinfection possible. | By suspect of infection Monitoring of IgG antibodies |

- negative, + positive

4. PERFORMANCE CHARACTERISTICS

Specificity / Sensitivity

100 samples IgA and 160 samples IgG and IgM were tested parallel in ELISA ANTI - VZV IgG/IgM /IgA and comparison methods (IFA/ELISA). The specificity and sensitivity are based on the results found.

Specificity: IgG 100% **Sensitivity:** IgG 100%
Specificity: IgM 99,4% **Sensitivity:** IgM 100%
Specificity: IgA 97,8% **Sensitivity:** IgA 66,6%

For the calculation of the sensitivity and specificity of the ELISA test, equivocal results were defined as positive results. The results refer to the groups of samples investigated.

Precision and reproducibility

Intra-assay reproducibility was determined by testing samples of different levels of antibody reactivity for at least 22 times in one test run. The coefficient of variation (CV) of the reactive IgG, IgM and IgA samples was < 10%.

Inter-assay reproducibility was determined by testing samples of different levels of antibody reactivity in 10 different test runs. The CV of the reactive IgG and IgM samples was < 10%.

Cross- reactivity

Cross-reactivity between antibodies to other Herpes viruses (HSV 1, HSV 2, CMV and EBV) cannot be excluded. Positive and equivocal results should therefore be confirmed with appropriate tests to exclude an infection with other herpes viruses.

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6. KIT COMPONENTS

1. **MICROTITERSTRIPS** [MTS]
One microtiterplate is supplied which contains 12 microtiterstrips of 8 breakapart wells. The wells are coated with purified, inactive antigen. Strips are colour-coded.
2. **PEROXIDASE CONJUGATE** [CONJ|POD]
One vial containing 12 ml of ready-to-use anti- human IgA, IgG or IgM Peroxidase Conjugate. Peroxidase conjugate contains 0.049% Thimerosal as preservative.
3. **NEGATIVE CONTROL** [CONTROL|-]
One vial of 1.2 ml containing human serum with 0,095% sodium azide as preservative. Ready to use.
4. **CUT OFF CONTROL** [CUTOFF]
One vial of 1.2 ml containing human serum with 0,095% sodium azide as preservative. Ready to use.
5. **POSITIVE CONTROL** [CONTROL|+]
One vial of 1.2 ml containing human serum with 0,095% sodium azide as preservative. Ready to use. The titre is reported on the label. Is no titre indicated on the label, the [CONTROL|+] is a positive reference sample and has a ratio $\geq 1,500$.
6. **TMB SUBSTRATE** [SUBS|TMB]
One vial containing 13 ml of ready-to-use tetra-methylbenzidine (TMB) substrate.
7. **SAMPLE DILUENT** [SPE|DIL]
One bottle containing 100 ml (2x50 ml) of ready-to-use sample diluent buffer. The buffer includes 0.049% Thimerosal as preservative.
8. **WASH SOLUTION 25X** [WASH|BUF|25x]
One bottle containing 80 ml (2x40 ml) of wash solution concentrate.
9. **STOP SOLUTION** [SOLN|STOP]
One bottle containing 15 ml of 0.95N H₂SO₄ stop solution. ready-to-use.

The safety data sheet (MSDS) is available upon request.

7. STORAGE AND STABILITY

Store all reagents at 2-8°C. Protect them from intense light and do not freeze. The expiration date of each component is indicated on the respective vial label. Do not use reagents beyond the expiration date. After opening, [MTS] must be stored at 2-8°C in the plastic bag with desiccant and are stable up to 4 weeks. The diluted [WASHBUF] is stable up to 4 weeks when stored at 2-8°C. Use only [MTS] with an intact vacuum packaging.

8. MATERIALS AND EQUIPMENT REQUIRED BUT NOT PROVIDED

- Test tubes for sample dilution
- Timer
- Micropipettes, multipipettes 10-1000 µl
- One-liter graduated cylinder, distilled water
- ELISA washer or multichannel pipette
- Spectrophotometer for micro-plates (450 nm/ reference wavelength 630/620 nm)
- Paper towels, pipette tips

9. WARNINGS OR PRECAUTIONS SAFETY PRECAUTIONS

The ELISA test is for [IVD] use only.

1. Only qualified and well-trained employees should carry out the assay procedure.
2. The instruction for use describes the applicable test method. In case of modification or applications others than the intended use, or the use of automatic processors, the user has to validate the procedure and take the responsibility for it.
3. Do not mix lot specific reagents, such as [MTS], controls and [CONJ|POD] from different kit lots. The [SUBS|TMB] doesn't have to be from the original test kit, but the lot of the [SUBS|TMB] has to be the same as indicated on the kit label. The [SPE|DIL] (except Immunocapture assays), the [WASHBUF|25x] and the [SOLN|STOP] can be used for all ELISA tests.
4. Seal all bottles properly after use in order to avoid bacterial contamination. All samples and kit components should be considered potentially infectious. All control samples have been tested for Hepatitis Bs antigen, anti-HIV I and II, anti-HCV (CE/FDA) and found to be negative.
5. The [MTS] are coated with inactive antigen. However, normal laboratory precautions should be maintained when handling with infectious material. Do not pipette by mouth.
6. Avoid contact with skin and mucous membranes when handling reagents, which contain preservatives (see kit contents). Wash thoroughly with water in case of contact and possibly look up a doctor.
7. Controls containing sodium azide may react with lead and copper plumbing, building up explosive metal acids. Flush with sufficient water when disposing of reagents.
8. The [SOLN|STOP] (0,95 N H₂SO₄) contains sodium hydroxide which may irritate skin and mucous membranes. Wash thoroughly with water in case of contact.
9. For disposal the legal regulations have to be followed.

10. SPECIMEN COLLECTION AND STORAGE

1. Microbially contaminated specimen may cause interference.
2. Lipaemic, hemolytic or icteric samples should only be tested with reservations although in our testing no negative influence has been found.
3. Suitable specimens are serum or plasma (heparinized, EDTA) samples obtained by standard laboratory techniques.
4. The samples should not be heat-inactivated since non-specific results may occur.
5. For ethical reasons, the use of CSF samples was only evaluated by testing pairs of CSF and serum samples (only if available), distributed by the WHO Collaborating Centre for Quality Assurance and Standardization in Laboratory Medicine in Germany. Therefore the use of CSF samples should be proved by testing a sufficient number of samples.
6. Patient samples should be stored at 2-8°C. For long term storage -20°C or lower is recommended. Avoid repeated freeze-thaw cycles.
7. **Note:** Diluted patient samples must be used on the same day.
8. **For the IgM specific antibody determination** an absorption of IgG / RF with IgG-RF-Sorbent must be performed in order to avoid interferences with high IgG titres and rheumatoid factors (RF).
(ALPHADIA IgG-RF-SORBENT [REF] VSB 100).

11. ASSAY PROCEDURE REAGENT PREPARATION

Bring all reagents to room temperature prior to use !

[WASHBUF] : Dilute the [WASHBUF|25x] 1:25 with distilled water
e.g. add 40ml of [WASHBUF|25x] to 960 ml distilled water and mix well.

Dilution of samples IgG / IgA: Dilute patient samples 1:101 with [SPE|DIL] e.g. 10µl sample + 1ml [SPE|DIL]; mix thoroughly.

For the **IgM specific antibody determination** the final dilution of 1:101 must be considered when an absorption with IgG-RF-Sorbent is performed. Please follow the instruction for the use of the IgG-RF-Absorbent used in the test.

Controls are ready to use.

Note: To take into consideration pipetting time, it is recommended that the [CUTOFF] is repeated after every 4 [MTS] (resp. after a pipetting time of >=5 min.) to evaluate the following patient's tests with the new calculated cut-off value. In case of a semiquantitative determination the [CONTROL+] should also be repeated the same way the [CUTOFF] was dispensed.

Take the required [MTS] out of the foil packets and place them in the holder. Possibly remaining wells of a [MTS] have to be stored at 2-8°C tightly sealed in the plastic bag provided, with the desiccant inside.

12. PIPETTING AND INCUBATION STEPS

- A. Pipette 100µl of the controls or diluted patient sample into the wells. Pipette 100 µl of sample diluent into well A1 (Blank).
- B. Incubate the wells at room temperature (21-25°C) for 30 minutes, protected from intense light.
- C. Wash the wells four times as described in section k. WASHING PROCEDURE
- D. Add 100µl of ready-to-use peroxidase conjugate to each well.
- E. Incubate the wells at room temperature (21-25°C) for 30 minutes, protected from intense light.
- F. Repeat washing as in section C above.
- G. Add 100µl of ready-to-use TMB substrate to each well.
- H. Incubate the wells at room temperature (21-25°C), in the dark for 10 minutes
- I. Add 100µl of stop solution to each well. Tap gently to ensure homogenous color distribution and read within 10 minutes.
- J. To read the plate, make sure that the bottom is free from moisture and that no air bubbles are in the wells. Read the absorbance of the well contents at 450nm on a suitable plate reader. On readers equipped with a dual wavelength facility set the reference filter to 620/630 nm.

Attention: The absorbance (OD) of the Blank must be always subtracted from the OD values of the controls and samples.

PROCEDURAL NOTES

Do not allow the wells to dry out between incubations.
Comply with the given incubation temperatures and times.

k. WASHING PROCEDURE

The washing procedure can be done manually with a multichannel pipette or on an automatic plate washer. Empty the wells, invert and tap dry on paper towel. Wash four times with a soaking time of approx. 30 seconds (300 µl).

13. SUGGESTIONS FOR TROUBLESHOOTING

In case that the ELISA instructions are followed strictly, the reagents are handled with care and the samples and reagents are pipetted carefully, the following kinds of errors can be avoided to a large extend.

| ERROR | POSSIBLE CAUSES |
|--|--|
| No colourimetric reaction after addition TMB substrate | No Peroxidase conjugate pipetted, contamination of Peroxidase conjugate (possibly with control sera during pipetting) may cause an inactivation. |
| Generally too high reaction | Incorrect Peroxidase conjugate (i.e. not from original test kit), incubation time too long or incubation temperature too high, water quality for Washing Solution insufficient (low grade of deionization) |
| Generally too weak reaction | Incorrect Peroxidase conjugate (i.e. not from original test kit), incubation time too short, incubation temperature too low |
| Reagent blank too high | Incorrect pipetting of sample diluent, contaminated reagents, reagents expired, exceeding of incubation time and temperature, external contamination of the bottom of microtiterstrips, (clean carefully!) |
| False positive / negative samples | Incorrect dilution of samples, microbially contaminated specimen |
| Unexplainable outliers | Contamination of pipettes, tips or containers or with metals (iron, copper etc.), insufficient washing |
| High variation (within a series) | Reagents (including microtiterstrips) not pre-warmed to room temperature prior to use. Washer is not washing correctly! |
| High variation (from series to series) | Incubation conditions not constant (time, temperature) high variation of incubation temperature, controls and samples are not carried out at same time (same intervals) check pipetting order, person related variation, strips dried out after washing (unreproducible results) |

14. VALIDITY OF THE ASSAY

All controls should be carried out with every test run.

The test must comply with the following validation criteria:

- OD-value of Negative Control should be < 0.100 ,
- OD-value of Cut-off Control should be > 0.200 ,
- Ratio of Positive Control/ Cut-off value should be ≥ 1.5
- OD-value of the Blank should not be higher than 0.100 .

If controls give invalid levels then results from the test samples are invalid too and retesting is required.

15. CALCULATION OF RESULTS

A QUALITATIVE CALCULATION

Calculation of [Cut-off Value]

The **Cut-off Value** is calculated from the absorbance of the Negative Control and the absorbance of the Cut-off Control and defines the Cut-off range.

$$\text{Cut-off Value} = \text{OD of the Negative Control} + \text{OD of the Cut-off Control}$$

$$\text{CUT-OFF RANGE} = \text{CUT-OFF VALUE} \pm 10\%$$

Interpretation of sample results:

| RESULT | DEFINITION |
|-------------------|--|
| negative - | OD value sample $<$ Cut-off value -10% |
| equivocal | OD value sample \geq Cut-off value -10% OD value sample \leq Cut-off value +10% |
| positive + | OD value sample $>$ Cut-off value +10% |

Equivocal results should be retested. Following the confirmation of the equivocal result the monitoring of the patient's antibodies is recommended in order to exclude unspecific reactions resp. cross-reactivity, which may also cause equivocal results.

B CALCULATION OF RATIO (CUT-OFF INDEX, COI):

Patient samples may also be quantified and interpreted by means of the calculation of the ratio (Cut-off Index, COI):

$$\text{COI} = \text{OD value of sample} / \text{Cut-off value},$$

whereby a ratio of 1.000 is equivalent to the Cut-off value.

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A semi-quantitative diagram is enclosed. The first point on the curve is obtained from the **Cut-off value** (y-axis) and the cut-off titre 1:100 (x-axis). The second point of the curve is obtained from the absorbance of the [CONTROL+] (y-axis) and their titre (x-axis) as indicated on the label. Drawing a straight line between the two points produces the semi-quantitative curve. The titre of the patient samples may be read from the curve. The graph is linear up to the titre of the [CONTROL+].

Samples with titres higher than the titre of the [CONTROL+], should be diluted further with [SPE]DIL according to the expected titre.
 For calculation of results, the dilution factor should be taken into consideration.

The calculated titres of the patient samples may also be (determined, shown, specified) as VU (VIRO-Units), e.g. a titre of 1:250 is equivalent to 250 VU.

D REFERENCE-UNITS

The **IgG** specific antibody activity of a patient sample may also be quantified by means of calculation of International Units. The quantitative test results (m IU/ml) are based on the **VZV Standard of the WHO** (WHO International Reference Preparations; Anti-varicella zoster immunoglobulin, human. 1 st. International Standard 1987). The Cut-off value correspond to **250 m IU/ml**. The [CONTROL+], i.e. the upper detection limit of the ELISA ANTI - VZV - IgG, is about **2000 m IU/ml**.
 [CUT_OFF] = **250 mIU (WHO)**

$$[\text{CONTROL}+] = \mathbf{2000\ mIU\ (WHO)}$$

If no titre indicated on the label, the [CONTROL+] is a positive reference sample and has a ratio $\geq 1,500$.

The [CONTROL+] for ANTI-VZV . IgA is a *positive reference sample* with a ratio ~ 1.500 . A

Quantification of the test result by means of titre specifications is therefore not possible.

Regarding diagnostic relevance and interpretation of results see **page 1**.

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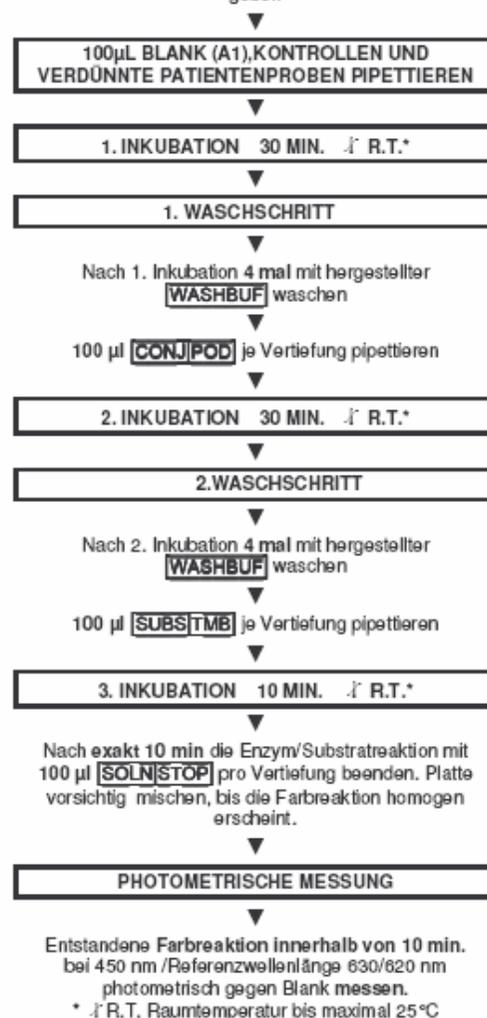
Prepare reagents


Dilution of samples:

Dilute patient samples 1:101
e.g. 10 µl sample + 1 ml **SPE|DIL**.

In case of pretreatment of the samples with RF Absorbent, the final dilution of 1:101 must be considered. Please follow the instruction for use of the RF Absorbent used in the test.

Controls/ calibrators are **RTU**!

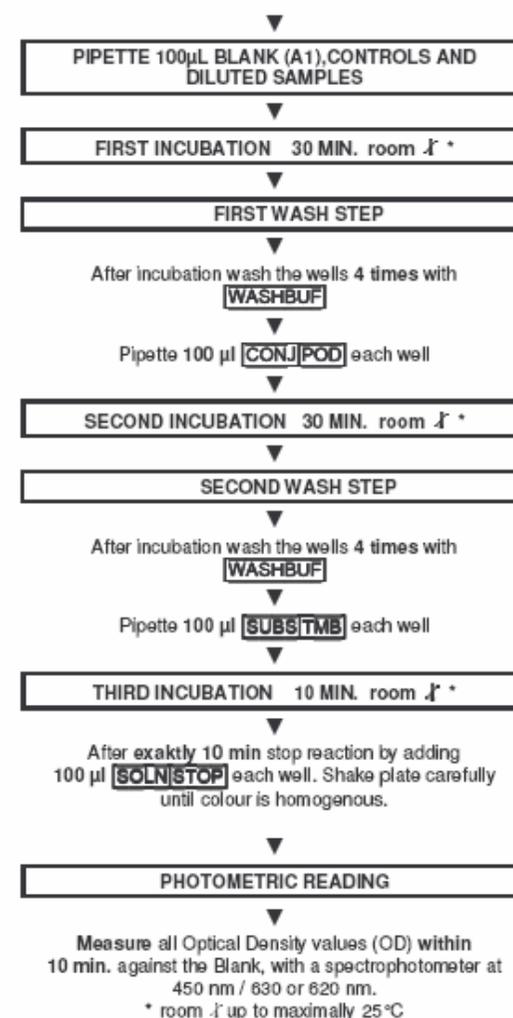Take required **MTS** out of the foil packet and place them in the holder. Do not touch the wells. Unrequired **MTS** should be put in the plastic bag with desiccant and closed carefully

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